


# Vascular Treatment Filter

using the Universal IPL with OPT™

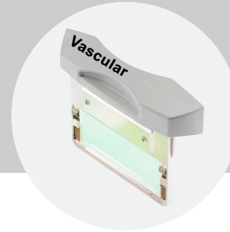

- **NEW** - Notch filter (530-650 & 900-1200 nm) especially designed for fine telangiectasia treatments
- **NEW** - Advanced Optimal Pulse Technology for precise tailoring of treatment parameters


## Clinical Results



Courtesy of Matteo Tretti Clementoni, MD

## The Lumenis OPT™ advantages

- Control pulse shape by equalizing energy distribution over the entire pulse or sequential pulses
- Gentler, more comfortable, patient-friendly procedures, with lower effective fluences
- High peak power, shorter pulses – optimal for IPL skin treatments using photorejuvenation and treatment of benign pigmented lesions

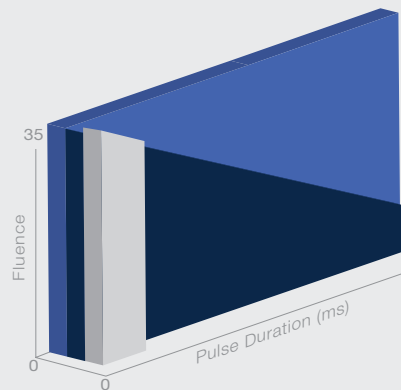

- M22 IPL with OPT
- High-powered IPL without OPT
- Limited capabilities of lower-powered systems

## Specifications using the Universal IPL

Wavelength	Vascular (530-650 & 900-1200 nm)
Fluence	10-48 J/cm <sup>2</sup>
Pulse Duration	4-20 ms
Pulse Delay	5-150 ms

Pulse characteristic	Multiple Sequential Pulsing
Repetition Rate	Up to 1 Hz
Spot Size	35 x 15 mm <sup>2</sup> ; 15 x 8 mm <sup>2</sup> ; 6 mm round
Cooling	Continuous contact cooling

© 2015 All Rights Reserved. The Lumenis Group of Companies. PB-2001373 Rev A. Lumenis, its logo, M22, OPT, Multiple-Sequential Pulsing, ExpertFilters, SapphireCool lightguides, are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice. Lumenis is ISO 13485:2012 certified. M22 is CE approved. | M22 is FDA cleared.

### AMERICAS

San Jose, CA, USA  
T +1 408 764 3000  
+1 877 586 3647  
F +1 408 764 3999

### EMEA

Dreieich Dreieichenhain, Germany  
T +49 6103 8335 0  
F +49 6103 8335 300  
Roma (RM), Italy  
T +39 06 90 75 230  
F +39 06 90 75 269  
Hertfordshire, UK  
T +44 20 8736 4110  
F +44 20 8736 4119

### JAPAN

Tokyo, Japan  
T +81 3 6743 8300  
F +81 3 6743 8301

### ASIA / PACIFIC

Beijing, China  
T + 86 10 5737 6677  
Gurgaon, India  
T + 91 124 422 07 95  
Kowloon, Hong Kong  
T + 852 217 428 00  
F + 852 272 251 51

For more information:  
[www.lumenis.com](http://www.lumenis.com)



[LUMENIS.COM/AESTHETIC](http://LUMENIS.COM/AESTHETIC)